




PRESSALIT®

Product Item number	Shower and changing table, electrically height adjustable R8592313021
Color Product description	<p> EAN:5708590357522 FDA:INQ/D124319 </p> <p>Lead grey</p> <p> Shower and changing table 3000, powered, height adjustable 27.6"; from 11.8" - 39.4", with wired hand control Foldable. Incl. safety rail and water collection tray with integrated outlet. Surface length 51.2", outer measurements 40.4" x 55.5". Excl. mounting kit. See Technical specifications (pdf) for information on mounting kits suitable for various wall types. <i>Replaces R8516 & R8534.</i> </p> 
Adjustability Adjustment height Product capacity Speed Noise level	<p>Height adjustable</p> <p>700 mm</p> <p> Maximum load: 440 lbs. When the table is loaded with 440 lbs, each mounting bolt supports 57 lbs. <i>This product has been tested statically with a load of 1100 lbs as the IEC 60601-standard requires testing with a safety factor.</i> </p> <p> Speed up: 15 mm/s Speed down: 20 mm/s </p> <p>48 dB</p>
Materials	Aluminium, polyurethane foam, plastics, powder-coated steel.

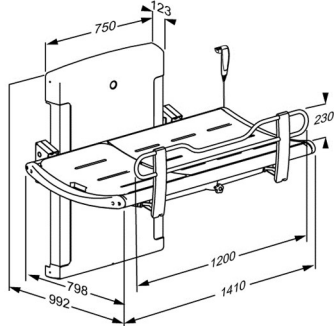

Product	Shower and changing table, electrically height adjustable
Item number	R8592313021
Ergonomy	<p>Pressalit offers a range of high quality, fixed or power adjusted changing tables designed to provide a stable platform suitable for changing adults and children in a variety of settings: schools, day care centers, institutions, hospitals and private homes.</p> <p>The showering and changing table 3000 series provides maximum safety and comfort. Besides an award winning design, the series includes comfortable and durable surfaces, including polyurethane foam and aluminum.</p> <p>The flexibility of the changing table provides optimal working conditions for the caregiver as well as a high level of comfort and safety for the user.</p>
Features	<p>A water collection tray with integrated water outlet prevents water from splashing on caregiver or the floor. The flexible hose attached to the center mounted outlet of the water collection tray, can be connected to either an outlet in the floor or on the wall. Adjustable neck support which can be moved to either end of the 51.2" changing tables.</p> <p>The slightly curved form makes the table extremely comfortable and safe.</p> <p>The safety rail can be folded easily up or down.</p> <p>The integrated pneumatic cylinders enable the table to be folded out of the way when not in use to maximize floor space.</p>
Mains connection (V)	120
Actuator (motor)	LA40CS Careline (Linak) 4000 N, IPX6
Control box	SMPS14 Careline (Linak), IPX6
Wired remote control	HB3 (Linak), IPX6
Power consumption (maximum load)	Power consumption by maximum load: 144 W Power consumption by standby: 2 W
Warranty	<p>Pressalit makes a commitment to remedy material and manufacturing defects for a period of three years from the invoicing date. NOTICE If the product is installed in an atmosphere heavy in salt or chloride, the warranty of the nursing bench will be reduced from the normal 3 years to 1 year. Salt and chloride have a corrosive effect on the products components and will cause a decrease in the life expectancy of the product. The damaging effects can, however, be reduced if the product is showered with tap water on a daily basis.</p>



PRESSALIT®






Product	Shower and changing table, electrically height adjustable
Item number	R8592313021
Cleaning	<p>This product should be cleaned with hot water or a mild cleaning solution. Do not use chlorine based, corrosive or abrasive cleaners as it may result in the surface being damaged.</p> <p>The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container.</p> <p>DISINFECTANTS: CaviCide</p> <p>The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.</p> <p>DECALCIFICATION OF ALUMINUM PARTS: Pressalit recommends that water deposits (Calcium, lime etc.) be removed from wall tracks and other aluminum parts using a solution of three parts of 30% acetic acid to seven parts tap water. When the aluminum parts have been washed down with this solution, the parts should be wiped over with a damp cloth wrung out in plain water.</p>








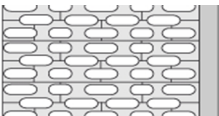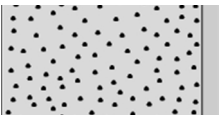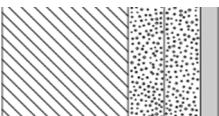



Product	Shower and changing table, electrically height adjustable
Item number	R8592313021

Mounting in different wall types

	Concrete >B25 V8608 - Mounting kit, 6 lag bolts Ø10 x 70 - plugs SX12. To be ordered separately.
	Solid Brick >= Mz12 V8608 - Mounting kit, 6 lag bolts Ø10 x 70 - plugs SX12. To be ordered separately.
	Solid Sand Lime Brick >= Hlz12 V8608 - Mounting kit, 6 lag bolts Ø10 x 70 - plugs SX12. To be ordered separately.
	Aerated Concrete P4 V8658 - Mounting kit, 6 lag bolts Ø10 x 70 - plugs GB14. To be ordered separately.
	Drywall (2 x 12,5 mm) water-resistant, re-enforced with wood V8608 - Mounting kit, 6 lag bolts Ø10 x 70 - plugs SX12. NB! Plugs are not needed. To be ordered separately.

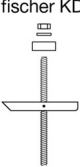

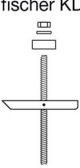PRESSALIT®

Product Item number	Shower and changing table, electrically height adjustable R8592313021
	<p>Drywall (2 x 12,5 mm) water-resistant, re-enforced with 22 mm plywood</p> <p>Wall fasteners (6 pcs) must be bought separately, for drywall it is possible to use e.g. SnapTight or WingIt anchors or similar.</p> <div>    </div>
	<p>Drywall (2 x 12,5 mm) water-resistant, re-enforced with steel plates (.08 inch)</p> <p>Wall fasteners (6 pcs) must be bought separately, for drywall it is possible to use e.g. SnapTight or WingIt anchors or similar.</p> <div>    </div>
	<p>Drywall (2 x 12,5 mm) water-resistant</p> <p>Mounting in this wall type is not recommended - risk of fatal injury to persons and walls.</p> 

It is the responsibility of the installer to ensure that suitable wall fasteners are used.



www.pressalit.com

Pressalit A/S
 Pressalitvej 1
 DK-8680 Ry
 Tel.: +45 8788 8788
 Fax: +45 8788 8789
 E-mail: pressalit@pressalit.com

PRESSALIT®