

| CAN-DAN ELECTRIC MOBILE HEIGHT ADJUSTABLE CHANGE TABLE Technical Specifications | |
|---|--|
| Measurements: | Table Sizes: 62x28, 72x28, 62x30, 72x30, 62x32, 72x32, 72x36, 72x40 Mattress: 2" thick Height Adjustment: 18" – 34" (floor to table top) Minimum Length for space: Length 62" table – 74", Length 72" table – 84" Lifting Capacity: 450 lbs (evenly distributed) Note: For the Can-Dan Mobile Changing Tables, the outlet should be no higher than 16" from the floor to the top of the outlet and preferably centred to the position of the table. |
| Product Group: | Change Tables |
| Description: | The Mobile Electric Height Adjustable Changing Table is a lifting device. Electrically powered hand control connected to a 24V actuator. The table is lowered for safe transfer for both user and caregiver. Mobile Electric Height Adjustable Change Tables are ETL Certified as a solution. Individual is secured using the 2 integrated safety belts. Height adjustable with wired hand control or optional foot pedal Suitable for use in schools, clinics, private homes, group homes, nursing homes, and other caregiving situations. |
| Materials: | Powder-coated steel frame. 5" locking casters. (4" locking casters on 36" and 40" tables Tabletop - Particle board with melamine. Safety belts – Woven fabric. 18 oz. vinyl mattress, mildew resistant, latex-free lead-free, phthalate free and fire retardant to NFPA 701 USA and CAN/ULC S109 standards. |
| Colour: | Frame – Grey Tabletop – White Mattress – Grey |
| Adjustability: | Height Adjustable. Optional fold down and swing away safety rails. |
| Ergonomics: | High quality powered change tables. Provides a stable platform for changing individuals. Ability to be lowered to a suitable level for safe transfer height. |
| Functionality: | Solid Surface with a 2" removable mattress. Level and straight to facilitate assisted or independent transfers. Hand control can be fastened to the front cover. Removable mattress fastened with Velcro. Caregivers use the hand control to lower the change table to a suitable height to safety transfer individual from their seating system to the padded tabletop. |
| Actuator: | LA31 Careline, (Linak) |
| Control Box: | CB9 Careline, (Linak) |




| Hand Control: | HB4 (Linak) |
|---------------|--|
| Motor: | Powered by 115VAC to 24 VDC power supply requires |
| | a GFI receptacle |
| Weight: | 300 lbs (depending on size) |
| Cleaning: | Approved disinfectants. Must be wiped down with a soft rag immediately after |
| | cleaning to prevent discolouration of the mattress. If not, the surface coating may |
| | deteriorate. |
| Warranty: | Can-Dan Rehatec Itd. as manufacturer warrants from the DATE OF DELIVERY TO THE FIRST CONSUMER that each item described below, as NEW AND UNUSED, will be free from any defects in material and/or workmanship for a period of: • 5 years on change Table Frame (Manufactured by Can-Dan Rehatec Ltd.) • 2 years on Linak Actuator • 1 year for accessories including: |

Picture & Diagram




Above: Change Table with optional Side-Rails



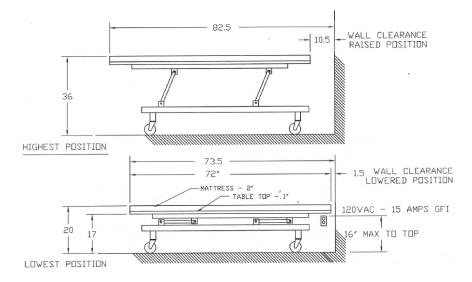


72" Mobile Specifications

62" Mobile Specifications

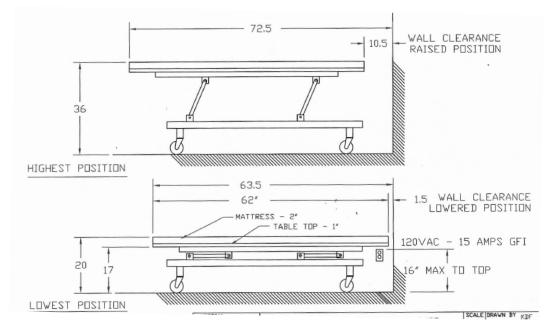



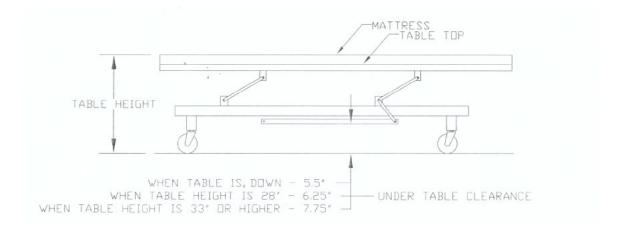

Under Table Clearance