

NEW
Early Intervention
Model Now Available



buddy
freedom to roam
roamer™

buddy
roamer™

The Buddy Roamer™ is a posterior walking aid that provides partial weight bearing postural support with mobility.

The ergonomically designed frame supports body weight while allowing the child to walk around freely and unaided.

more than
just a walker



Our Aims and Objectives

The opportunity to walk is incredible for our client group, but if, at the same time, we also see benefits of communication and relationships, we also affect that much more elusive thing we call 'quality of life'!

Encouraging Communication

As we set the unit up for Sarah-Jane her mother was asking 'what video do you want to see?' It was Sarah-Jane's first time to walk unaided and as she took her first tentative steps her Mum called out 'Sarah-Jane, come over here' but her daughter didn't listen. Instead she made her way across the room and very deliberately lifted, from the middle of the pile, the video she wanted to watch and handed it to her Mum!

This lovely, non-verbal 9 year old girl had made her wishes known - loud and clear! - Mobility allowed her to communicate!

Enabling Relationships

Imagine the excitement when a Dad said to us recently "This is the first time in ten years that my son was able to give me a hug!"

Experiencing New Possibilities

If you have ever seen the look on a child's face as it takes its first unaided steps, there is a magic that cannot be described! The Buddy Roamer™ has the potential of opening up new worlds not yet experienced!

Completely open to the front

There is no obstruction to interaction with other people or objects. This allows the child to have hands-free interaction with the surroundings to learn, explore and discover.

Ergonomic

The frame has been developed so as to allow easy access through doors and down hallways while still providing excellent stability for exploration both indoors and out.

"The Buddy Roamer™ fulfils my dream of a functional unit that enables children to walk following the natural gait pattern, and the fact that it is open to the front encourages communication, relationships, interaction, play and learning."

Girvan MacCorkell,
CEO of Moorings Mediquip



For more testimonials, inspirational stories, videos, features and much more visit www.buddyroamer.com

Now Introducing – the Early Intervention model – the Buddy Walkabout.



Mulholland
Positioning Systems Inc.




Lighter and smaller frame suits younger users.

Building on the heritage of 40 years of walkers, the Buddy Walkabout is today's ideal solution to providing walking assistance as soon as it is needed, with comfort and stability.

Free Lateral Movement

The Centre of Gravity of the Buddy Roamer™, being to the rear, means that there is little or no weight on the front castors allowing for effortless free movement and left to right weight shift during the stride phase.


This design feature means that at all times stability is maintained at the pelvis.

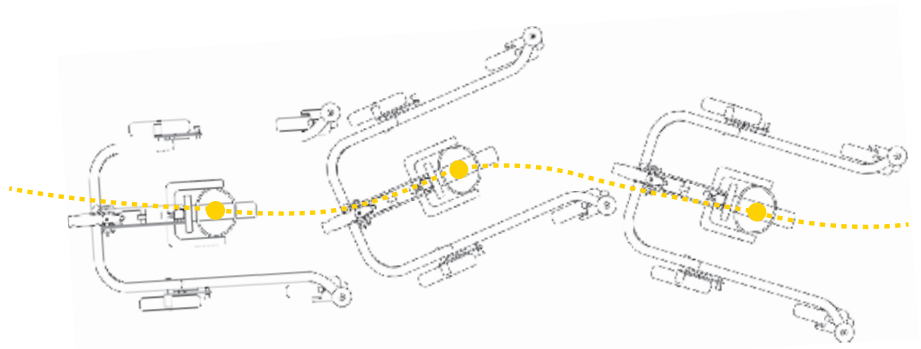

Adjustable Dynamic Spring

Not only does the Dynamic Adjustable Spring allow the user to naturally follow the gait pattern, it provides spring assisted lift during the stride phase. It should be noted that the angle of the spring also gives a forward component to the assistance.

The spring is easily adjustable to suit individual users, even when the user is in the unit.


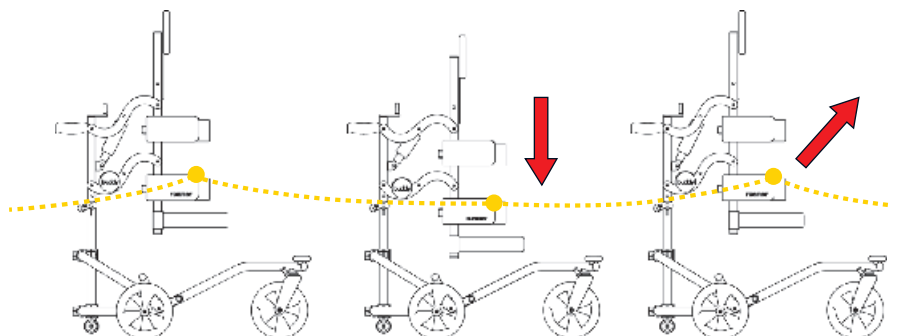


Product Sizing Range

| Measurements | Walkabout | | Size 1 | | Size 2 | | Size 3 | | Size 4 | |
|-----------------|-----------|---------|----------|---------|----------|---------|-----------|---------|-----------|---------|
| | cm | inches | cm | inches | cm | inches | cm | inches | cm | inches |
| Headrest Height | 52 – 112 | 20 – 44 | 70 – 130 | 27 – 51 | 80 – 140 | 31 – 55 | 110 – 200 | 43 – 78 | 110 – 200 | 43 – 78 |
| Inside Leg | 15 – 45 | 6 – 18 | 20 – 70 | 8 – 27 | 30 – 80 | 12 – 32 | 40 – 100 | 16 – 39 | 40 – 100 | 16 – 39 |
| Chest Width | 15 – 22 | 6 – 8.5 | 15 – 23 | 6 – 9 | 21 – 36 | 8 – 14 | 25 – 40 | 10 – 16 | 30 – 50 | 12 – 20 |
| Hip Width | 15 – 22 | 6 – 8.5 | 15 – 23 | 6 – 9 | 21 – 36 | 8 – 14 | 25 – 40 | 10 – 16 | 30 – 50 | 12 – 20 |
| Max User weight | 30kg | 66lbs | 50kg | 110lbs | 50kg | 110lbs | 100kg | 220lbs | 100kg | 220lbs |

For early intervention – custom supports available.

Height Range Available

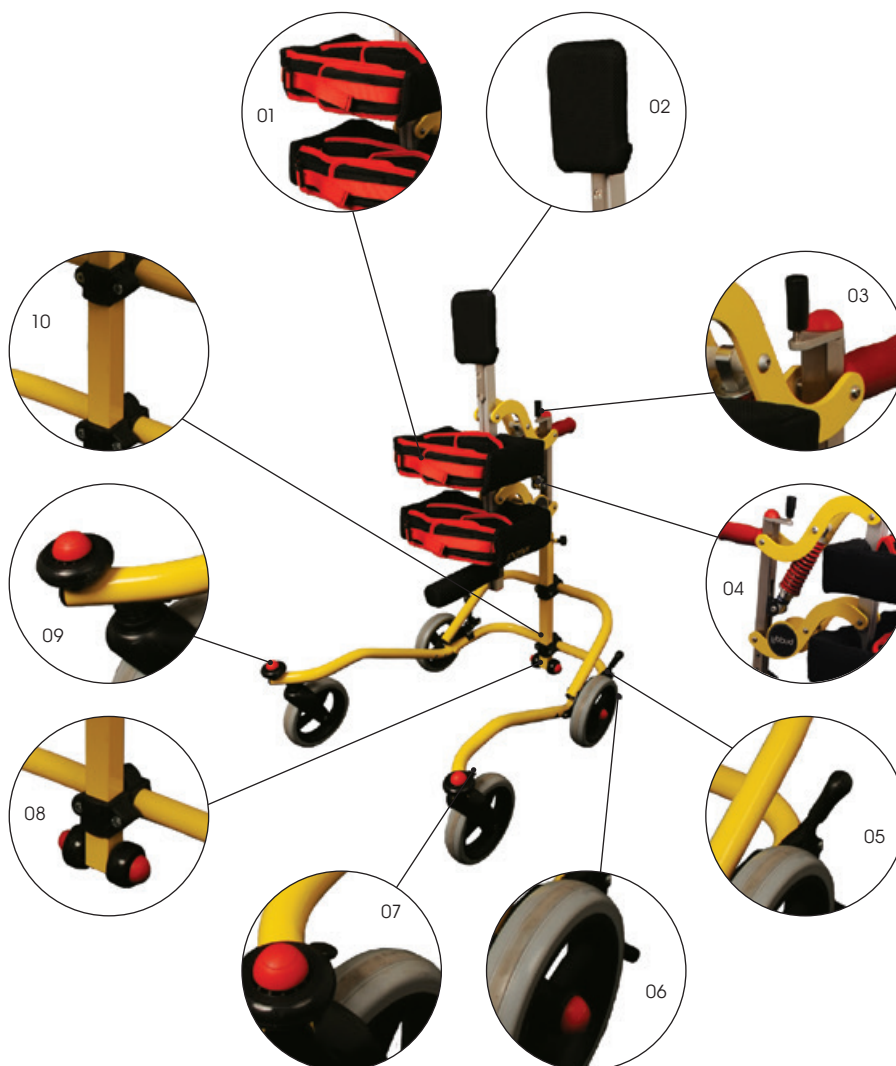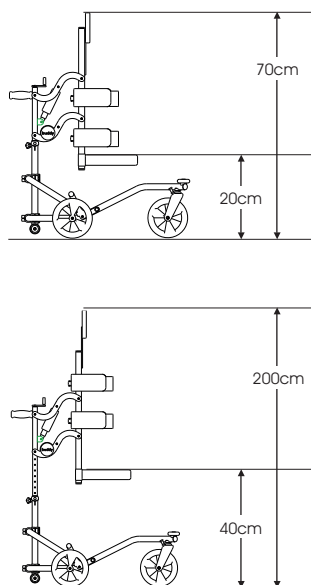

Features of a Buddy Roamer™

01. Adjustable Lateral Supports
02. Telescopic Headrest
03. Fine Height Adjustment
04. Adjustable Dynamic Spring
05. Parking Brake
06. Anti Roll-Back Lock (optional)
07. Directional Lock
08. Anti-Tip Wheels
09. Bumper Wheels
10. Prone Angle Adjustment

**buddy
roamer™**
www.buddyroamer.com



Moorings Mediquip
51 Slaght Road
Ballymena
Co. Antrim, BT42 2JH

Tel +44 (0) 28 2563 2777
Fax +44 (0) 28 2563 2272
www.mooringsmediquip.com

Authorised Dealer