



PROTAC®
INTEGRATING SENSES



Protac® SenSit® Sensory Learning Station

The Protac® SenSit® Sensory Learning Station is ideal for individuals with mental or motor agitation. It is designed to provide supportive seating and sensory input. The SenSit® is filled with balls in the seat and backrest, as well as in the side wings and neck cushions. The balls stimulate the proprioceptive sense and deep touch receptors through dynamic pressure, helping individuals feel secure, relaxed, and concentrated.

The Protac® SenSit® Sensory Learning Station is available in two sizes and two fabric options.

- Oxford is a flame-retardant, 100% polyester fabric for indoor use.
- Nature is a waterproof, wipeable, PU-coated polyester fabric for indoor and outdoor use, ideal where hygiene and protection from spills are a priority.



Protac® SenSit®

Suitable for all ages.

Oxford Fabric Colours:

Aqua · Dark Grey · Dusty Green · Lime · Off White



Nature Fabric Colours:

Navy · Khaki




Protac® SenSit® Junior

Recommended for users approximately 3-10 years old.

Oxford Fabric Colours:

Aqua · Dark Grey

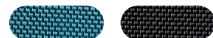

Nature Fabric Colours:

Navy · Khaki



Size	Floor to seat height	Floor to back height	Puff height	Base width	Weight capacity	Client age range
SenSit®	44 cm (17")	85 cm (33")	40 cm (16")	105 cm (42")	200 kg	All ages
SenSit® Junior	34 cm (13")	65 cm (26")	32 cm (13")	81 cm (32")	90 kg	Approx. ages 3-10

Protac® SenSit® Nature

The Protac® SenSit® Nature fabric is designed for indoor and outdoor use. Its durable, waterproof fabric, taped seams and zippers make it a practical choice for playgrounds, schoolyards and other outdoor spaces, providing familiar sensory support wherever it is used.

Indoors, the wipeable surface is ideal for environments with strict hygiene requirements and for users with incontinence. The easy-to-clean material helps protect against spills and moisture, making everyday care simpler.





PROTAC®
INTEGRATING SENSES

Key Features

Adjustable sensory support

The four sensory-stimulating wings can be positioned around the body to create a calm, enveloping feeling.

Dynamic pressure

Balls in the backrest, wings, and seat provide deep, dynamic pressure that stimulates the proprioceptive sense and sense of touch. This can promote body awareness and relaxation.

Adaptable seating

The ball-and-granulate filling makes the chair adaptable to each individual user.

Postural support

The wings can be positioned to help stabilize sitting posture while giving the user flexibility to adjust the level and location of sensory input.



V1.2.09/26

Accessories

Protac® SenSit® Puff

The SenSit® Puff can be used as a footrest or an additional seat. It has balls in the top and granulate filling in the base, providing additional sensory input and support. The SenSit® Puff also gives therapists, teachers, and caregivers a comfortable way to sit alongside the user at eye level, creating a setting for interaction and conversation.

The standard SenSit® configuration includes the SenSit® Puff.

Protac® SenSit® Frame

The Protac® SenSit® Frame makes the SenSit® chair easier to use for people who need extra support getting in and out. Handles assist with sitting down and standing up, while the forward tilt supports transfers and comfortable positioning. Rear wheels make it easier to move the chair between spaces.

The SenSit® Frame is not compatible with the SenSit® Junior.



Incontinence Liner for Oxford Fabric

An optional inner liner provides a protective barrier to the seat portion of the SenSit®. Available for both SenSit® and SenSit® Junior with Oxford fabric.